

KERATOR Attachments

For Root

There are two types for root retained overdenture.


The Post type is used in direct technique and there are straight, 10 degree, 20 degree angled females. The female part is made of Titanium alloy (TI6AL4V) and is coated with a Titanium Nitride coating.

And the Cast type has a pattern that is incorporated into wax-pattern for post coping.


The Root Kit includes 4 attachments (Male&Female), 3 kinds of burs (initial drill, main pilot drill, countersink diamond bur). The Diamond bur is used for the shoulder on the top of the root accept the KERATOR Female.



A. Direct Technique

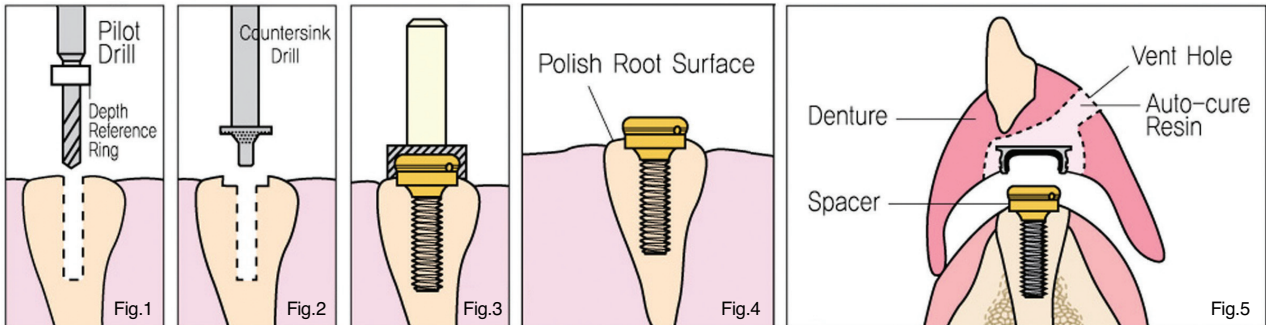

NON-DIVERGENT ROOT

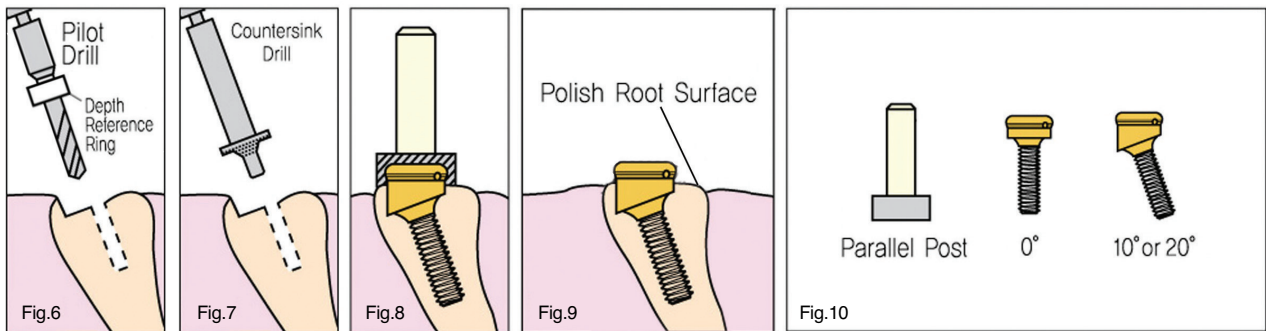

DIVERGENT ROOT

B. Cast-to Technique

