




ROTO QUICK

The versatile
system coupling
between hose
and instrument

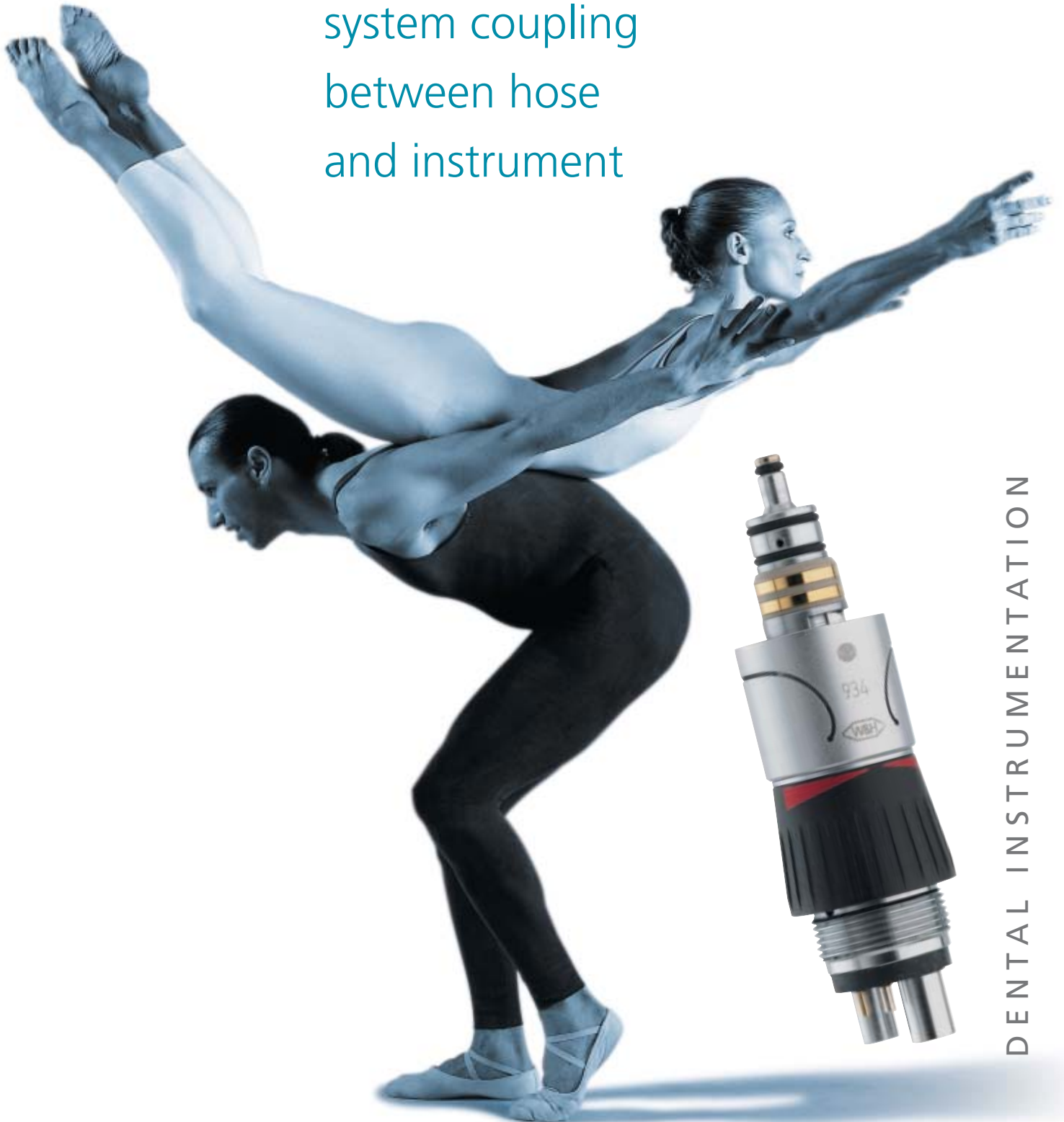

DENTAL INSTRUMENTATION

Roto Quick – one for all!

Fast and easy changing of instruments
is all part of a smooth treatment
session.

The Roto Quick system coupling is a versatile
interface between turbine hose and all
W&H instruments with the Roto Quick connection.
One quick
manipulation is all
it takes – and the required
instrument is ready to use!

The advantages of the Roto Quick are obvious:

- one system coupling compatible with most
W&H air instruments
- free and easy rotation
of instruments
through 360°
- fast and safe instrument
interchange
- non-retraction valve for optimum hygiene



One click – and ready to use!

The Roto Quick system coupling from W&H allows speedy instrument changing and thus continuous treatment.

It can be fitted on all standard turbine hoses. All W&H turbines, air motors scalers and polymerization lamps with Roto Quick connection will fit Roto Quick.

Important detail. The non-retraction valve prevents the back flow of contaminated spray water thereby minimizing the risk of infection.

The spray water flow can be smoothly adjusted with the control on Types 934, 914 and 915.



Smooth rotation. When fitted, instruments can be rotated effortlessly through 360° – with W&H instruments, this free rotation is especially noticeable when the instrument is in operation since the ring is mounted on ball bearings. No twisted hose, no drag on the wrist! This facilitates precise manipulation and contributes to fatigue-free working.

The instrument is conveniently removed by means of pressure points. W&H places great emphasis on operating convenience.



Guaranteed. The Roto Quick system coupling complies with the strict European product directives for medical equipment. It comes with a 1 year warranty.

ROTO QUICK

Technical data:



Roto Quick **with light**

Roto Quick **without light**



924 with 6 hole connection, ISO 9168



903 with 2 (3) hole connection, ISO 9168



904 with standard 4 hole connection, ISO 9168




934 with 6 hole connection, spray regulation, ISO 9168



914 with standard 4 hole connection, spray regulation, ISO 9168

915 with Siemens 4 hole connection, spray regulation

Roto Quick system: versatile and efficient

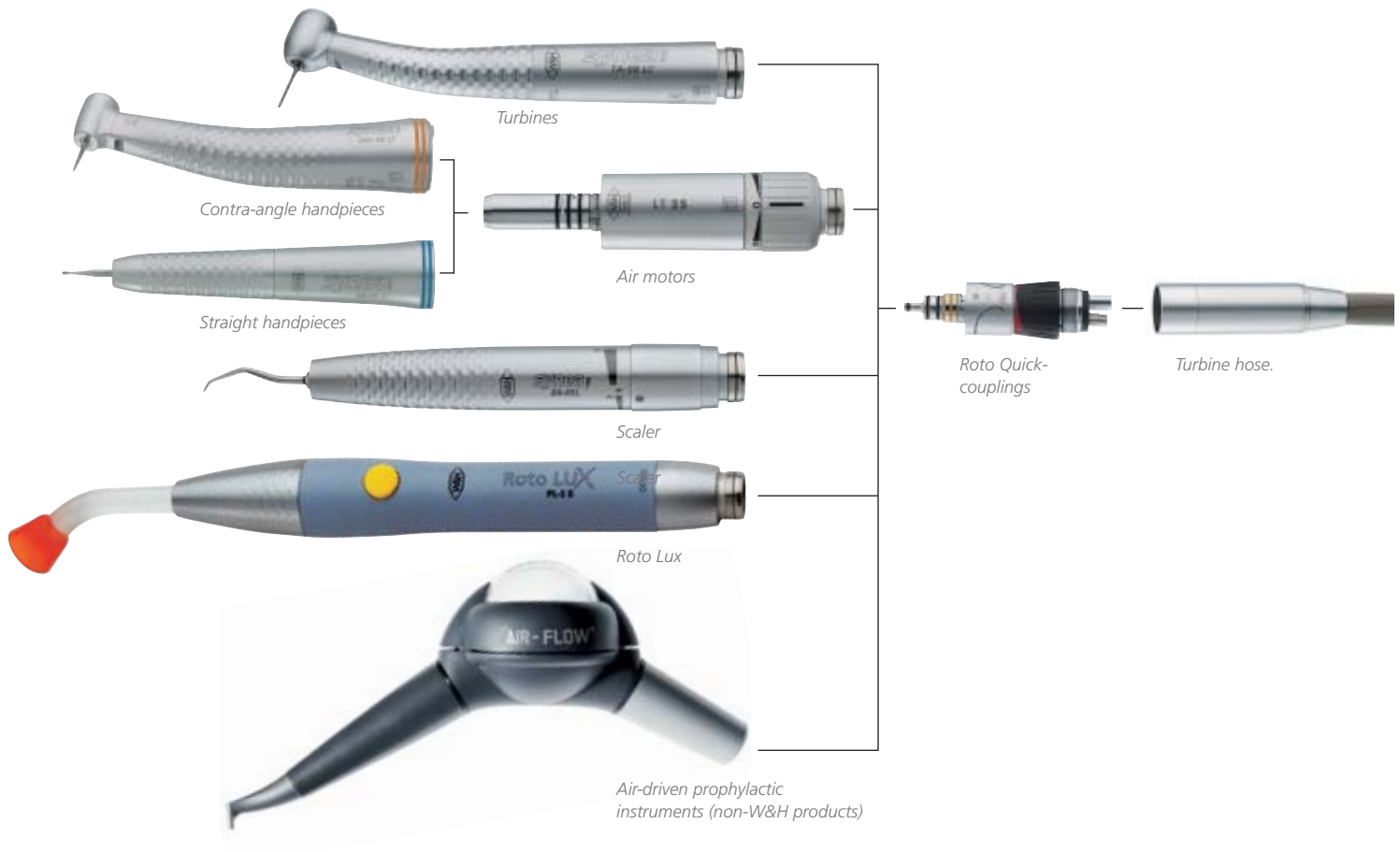

Technically perfect solutions from W&H.

Order W&H Roto Quick system couplings from your dealer.



DENTALWERK

W&H Dentalwerk Bürmoos GmbH
Austria, A-5111 Bürmoos

✉ Postfach 1, Ignaz-Glaser-Straße 53
☎ +43/62 74/62 36-0 @ office@wh.com
☎ +43/62 74/62 36-55 🌐 www.wh.com



ISO 13485
Annex II
93/42/EWG
EC 0297