

With Valmy, your confidence is also highly protected

As a resolutely high-tech company, Valmy chooses to have controlled industrial organisation. It uses European supply sources for their reliability, a French manufacturing site for proximity, machines designed to specifications for productivity, and continual control by an in-house laboratory to ensure total quality that is validated by an EU notified body. For our distributor partners and for all users of Valmy products, this guarantees complete control of the processes and costs, thereby ensuring that the range of products is competitive.



Valmy creates a new future for personal protection

As a 21st century company, Valmy observes the needs of the users and develops the most appropriate solutions. By gaining precise knowledge of practises we are able to ensure that each of our products is designed to satisfy the most demanding users' expectations. Our products are made a reality through the involvement and mobilisation of a multidisciplinary project team. They devise a comprehensive response that combines research, development, purchasing, production, marketing, publicity, advertising, sales and logistics. This approach complies with currently applicable standards, as part of our ISO 9001 accreditation.



Your confidence under high protection



**FFP2 respiratory protection mask
NHS Product Code BTP058**





FFP2 Respiratory protection masks

NHS Product Code : BTP058

Valmy® respiratory protection masks protect the user against the risk of inhaling infectious particles - such as dusts, pollution, virus, allergens...

Applications

| Maximum values permitted to comply with EN 149 : 2001 | | | | | | |
|---|--------------------------|----------------------|---------------------------------|----------|---------------------------------|--|
| Class | Filter media penetration | Total inward leakage | Inhalation breathing resistance | | Exhalation breathing resistance | Applications |
| | | | 30 l/min | 95 l/min | 160 l/min | |
| FFP2 | < 6% | < 8% | < 0,7 | < 2,4 | <3,0 | Antivirus (Avian flu, SARS, Tuberculosis), anti-pollution, anti-allergens. Dusts such as concrete, cement, plaster (demolition, carcass work, sanding), softwoods (sanding, cutting) fine particles of paints and resins, plastics (blasting, sanding) |

Features

- Latex free
- Hypoallergenic & soft
- Extra light : less than 6 gr
- Perfect hygiene : individually packed
- Adjustable nose piece
- Single use
- White non woven textile
- Two elastics for better fit
- Shelflife : 4 years from manufacture date
- Fluid resistant to 120 mmHg
- Bacterial filtration efficiency (Nelson) : > 99 %

Standards


Our Valmy masks meet the safety requirements laid down by the Council Directive 89/686/EEC. They fully meet all relevant European Standards EN 149 : 2001 by a notified body. Classification : FFP2 (FFP1 and FFP3 also available)

Caution

The user is responsible for the choice of the mask and must check that the mask is suitable (after risk assessment. Before use, please fully read and understand the user manual scrupulously.

Packaging

Individual packing

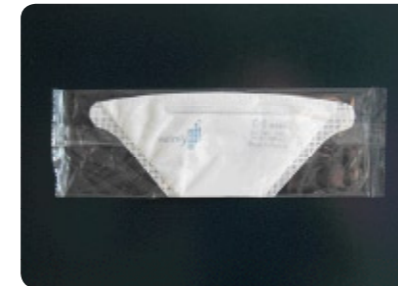

- Protection and hygiene guaranteed
- Marking according to regulations on the mask

Dispenser box



- Precise and fast identification of protection level reference, Lot no., shelf life...
- User information : applications, fitting instructions

Composition

- Filter media : polypropylene
- Elastics : polypropylene, polyamide, elastane
- Nose piece : steel & polypropylene
- Individual pack : polypropylene

Advantages

Exclusive design for :

- Greater filtration surface
- Leak-tightness increased
- Fits all faces
- Extra light - less than 6 gr

Limitations

This mask will protect you against dust, mists and aerosols in concentrations not exceeding 12 x MAC (Maximum Admissible Concentration).

Fitting instructions



1



2



3



4

ISO 9001 quality management system
Manufactured in France
More than 80 % of the mask is recyclable
100 % of manufacturing waste is recycled



Deltalyo & Valmy

Consultant manufacturer

VALMY
 Espace Valmy - Rue de Charlieu
 42300 ROANNE - FRANCE
 Tel. +33 (0)4 77 23 11 21
 Fax +33 (0)4 77 72 54 57
 info@valmy.eu

France - Export Sales

DELTALYO & VALMY
 105 Avenue Jean Jaurès
 69921 OULLINS CEDEX - FRANCE
 Tel. +33 (0)4 78 57 80 00
 Fax +33 (0)4 78 57 25 02
 info@deltalyo.com

UK Sales

DELTALYO & VALMY UK
 Units 21 & 23 Scott Road, Off Park Avenue
 Luton - Beds - LU3 3BF
 Tel. + 44 (0)1582 563600
 Fax +44 (0)1582 563323
 deltalyo.kestrel@deltalyo.com

www.valmy.eu