

4 reasons to use Luxator® RootPicker

#1

Rigid and slim

Designed to offer a strong grip with a slim design that is easy to insert into the alveolus.

#2

Unique design

Ideal for handling components used in surgery, such as implant accessories.



#3

Four corner grip


Provides an excellent grip on root tips.

#4

Straight or Angled

For best access.

Technique

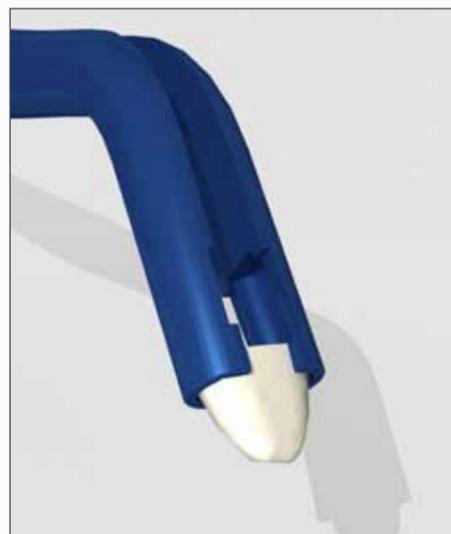

Luxator® RootPicker grips debris after extraction



During extractions annoying tooth fragments or loose root tips remain attached to the periodontal ligament. These pieces are very small and slippery, making it difficult to grip them with ordinary tweezers or forceps. RootPicker instruments address these problems with the necessary properties to complete the extraction successfully.

Luxator RootPicker is a short and rigid instrument that makes it easy to insert into the alveolus. The cleverly designed tip provides a sure grip on the debris.

Their straight and angled tips make them ideal for use in the anterior area as well as in the posterior maxillary and mandibular areas.

Luxator® RootPicker

Prod No	Product Description	Quantity
506450	Straight RootPicker	1 pc
506451	Angled RootPicker	1 pc
Kit		
506452	Kit K2 1 Straight + 1 Angled RootPicker	2 pcs

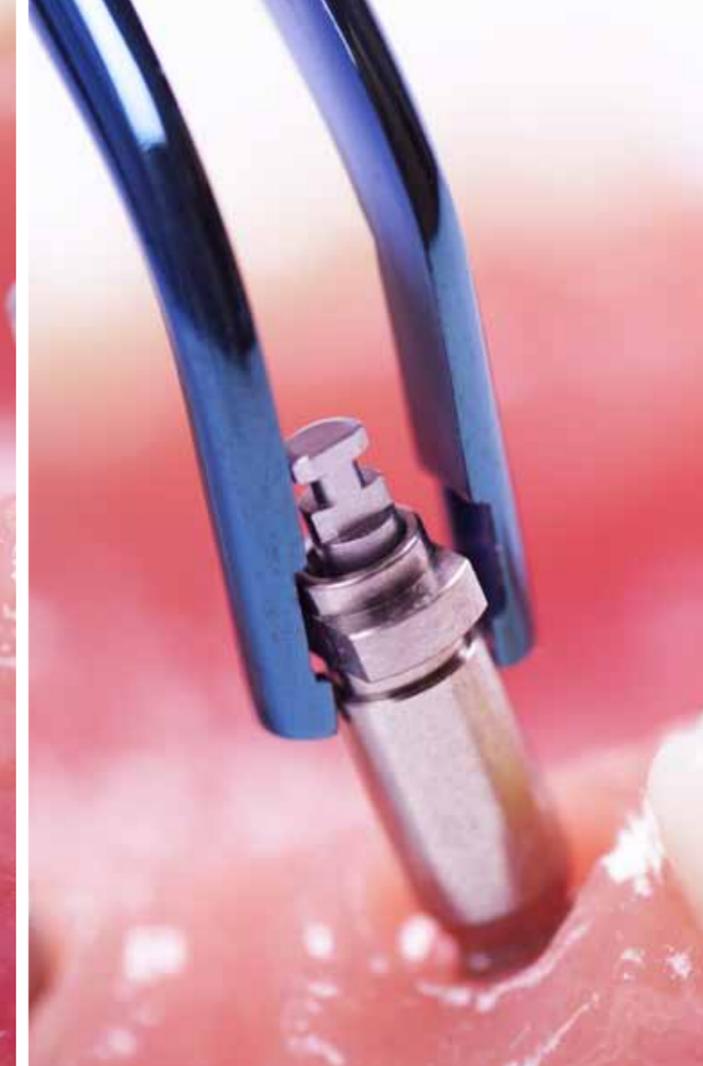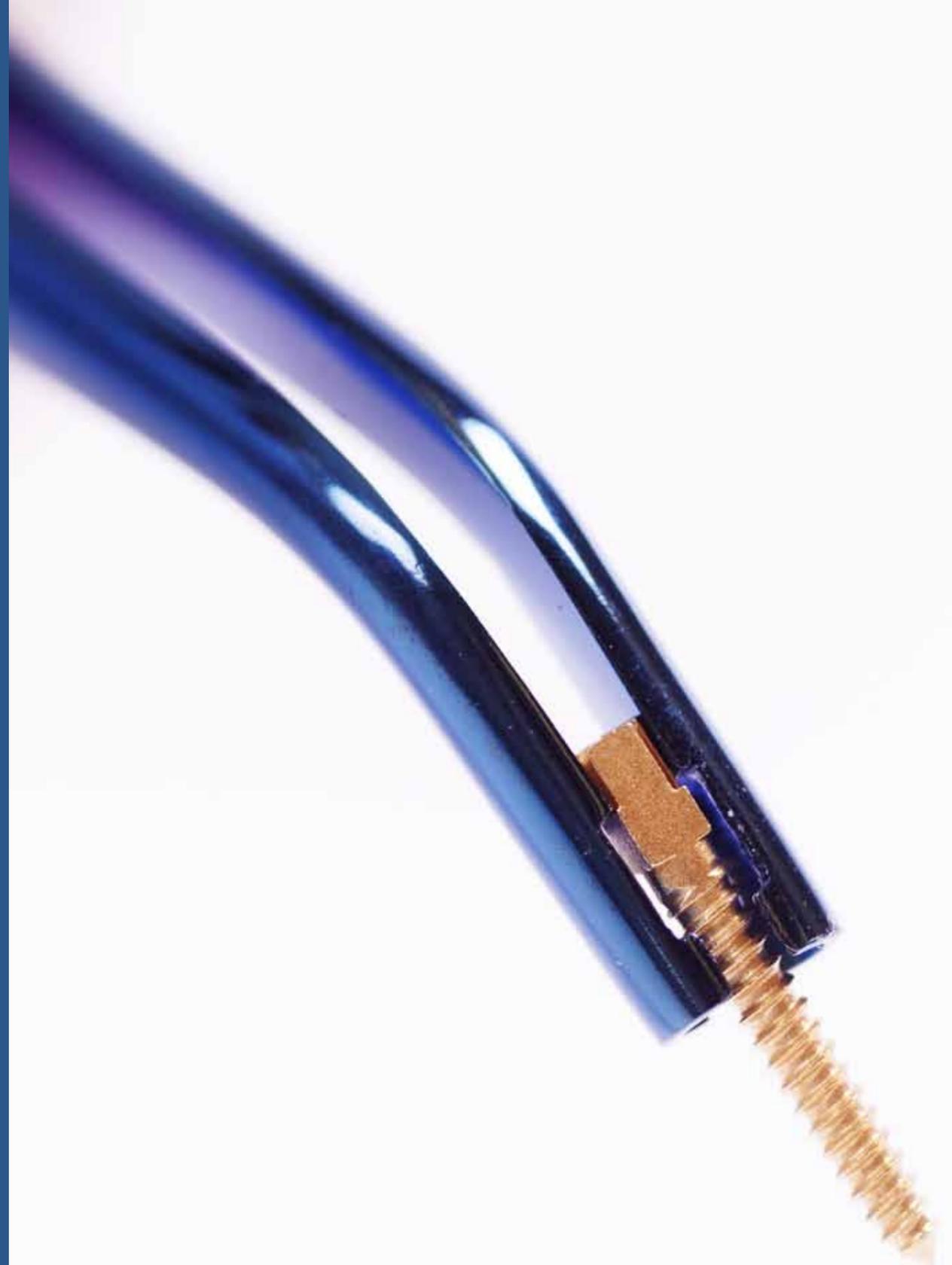

Luxator® RootPicker grips implantology components

RootPicker also provides excellent grip for the various components used in surgery such as healing abutments, parallelism indicators, screwdrivers and posts for immediate prosthesis.

Their straight and angled tips make them ideal for use in the anterior area as well as in the posterior maxillary and mandibular areas.

Luxator® RootPicker

Prod No	Product Description	Quantity
506450	Straight RootPicker	1 pc
506451	Angled RootPicker	1 pc
Kit		
506452	Kit K2 1 straight+1 angled RootPicker	2 pcs